Triad[®] DuaLine[®]

Visible Light Cure
Dual-Cure Reline Material

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MODE D'EMPLOI

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Indications for Use:

TRIAD DualLine Dual-Cure Reline Material is indicated for full or partial denture relines. TRIAD DualLine is a methyl methacrylate monomer-free material that combines both self-cure and visible light cure. TRIAD DualLine is a dual paste catalyst/base reline system dispensed from a cartridge.

Contraindications:

1. TRIAD DualLine is contraindicated for patients and users with a history of allergic reaction to a) urethane resins; and/or b) methyl methacrylate monomer (contained in TRIAD VLC Bonding Agent only).
2. TRIAD DualLine is contraindicated for patients with infected, inflamed or lacerated oral tissue. These conditions should be alleviated before treatment.
3. TRIAD DualLine is contraindicated for patients with excessive undercuts.

Warnings:

1. TRIAD DualLine and TRIAD VLC Bonding Agent contain polymerizable monomers which may cause skin sensitization (allergic contact dermatitis) or other allergic reactions in susceptible persons. Wash thoroughly with soap and water after contact. If skin sensitization occurs, discontinue use. If dermatitis or other symptoms persist, seek medical advice.

Precautions:

1. Users with special skin problems, cuts, or abrasions should wear protective gloves.

2. When grinding prosthodontic resins, proper ventilation, masks and vacuum systems should be used.
3. All cure times cited are appropriate for TRIAD 2000 and TRIAD II curing units. Consult your Technique Manual or call Technical Support for TRIAD I cure times.
4. Store TRIAD Materials at or below 24°C/75°F and away from direct sunlight. TRIAD Materials may be refrigerated to extend shelf life. **Always bring DualLine to room temperature before use.**
5. For best results, remove the denture 1½-2 hours prior to the reline procedure to allow sufficient tissue relaxation.
6. To avoid compromising occlusal relationships, prepare maxillary and mandibular relines sequentially.
7. If patient has inadequate salivary flow, tissue lubrication may be needed.
8. Always cool the denture to room temperature before applying DualLine. A warm denture will reduce working time and may cause excessive heat in the patient's mouth.
9. Extrude DualLine slowly to avoid excessive pressure build-up and cartridge leakage.
10. DualLine must undergo final cure in the TRIAD VLC unit to achieve full physical properties.
11. Thermoplastic denture bases relined with DualLine may be distorted in the TRIAD Unit.

12. Unused TRIAD Materials should be polymerized prior to disposal.

Adverse Reactions:

1. Allergic contact dermatitis and other allergic reactions may occur in susceptible individuals.
2. Dust will be generated when grinding TRIAD Materials. Eye, skin and respiratory irritation may occur if appropriate engineering controls are not used.
3. Some patients may experience an unpleasant taste from uncured TRIAD Materials.

STEP-BY-STEP INSTRUCTIONS

Check Denture Relationships

1. Before denture preparation, check and record the patient's vertical dimension of occlusion (VDO), and the anterior-posterior (A-P) denture position. A bite registration may aid in this process.

Denture Preparation

1. Reduce peripheral borders of denture 1 mm and contour as shown in Figure 1. TRIAD VLC Materials may be used to extend borders if needed.

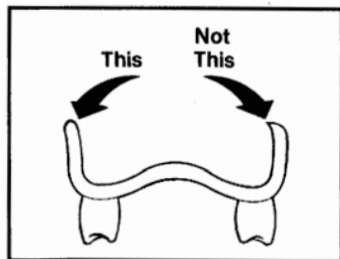
Figure 1 - Reducing/Contouring Peripheral Borders

2. **“Tissue Stops” are strongly recommended** to control the vertical dimension and thickness of the reline. See Figure 2.

NOTE:

Tissue stops are also recommended if the VDO is to be opened slightly. Roughen the denture surface to promote bonding and proceed with the technique establishing the new VDO.

Removable partial dentures should not require tissue stops. (See Reline Procedure, step 4.)

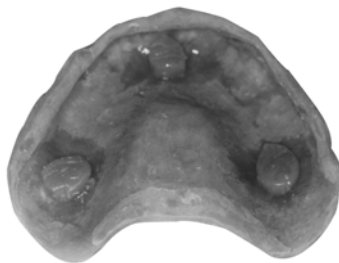
Figure 2 - Tissue Stops

- a. Reduce the denture surface by 1 - 2 mm.
- b. Apply TRIAD Bonding Agent to areas where tissue stops will be placed.
- c. *Bench set the denture for 2 minutes and cure in the TRIAD unit for 1 minute.*
- d. To form tissue stops, place small pieces of TRIAD Hi-Flow or Regular Reline Material in a tripod arrangement. (See Figures 3 and 4.) An alternative arrangement on a maxillary appliance may include

one tissue stop and a postdam palatal seal. (See Figure 5.)

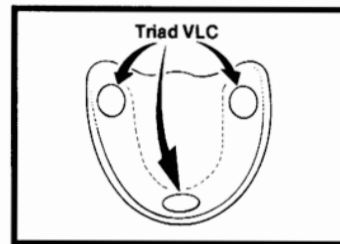
Figure 3 - Maxillary with Tissue Stops

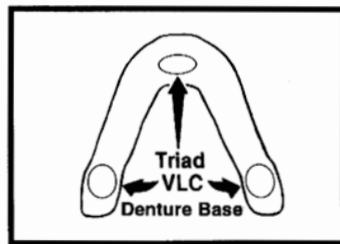
Figure 4 - Mandibular with Tissue Stops

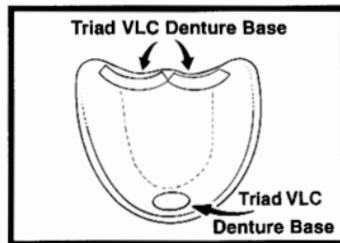
Figure 5 - Maxillary with Tissue Stop and Postdam

- e. Insert the denture. With light pressure, re-establish occlusion, VDO, and A-P relationships.
- f. Remove the denture and verify that the thickness of the tissue stops is 1-2 mm.

Remove any excess material and/or readapt if necessary.

- g. *Cure the tissue stops for 2 minutes in the TRIAD unit.*
3. Apply TRIAD Bonding Agent to all prepared areas including tissue stops. Bench set for 2 minutes and *cure in the TRIAD unit for 1 minute.*

The denture will be slightly warm from the tissue stop curing cycles. **Allow the denture to cool before proceeding.** The dental unit air syringe will quickly cool the denture.
4. Brush a thin coating of TRIAD Model Release Agent (MRA) on the denture teeth and buccal surfaces of the denture to simplify clean up.

RELINING PROCEDURE

NOTE:


Excessive pressure on the grip handle of the Dualine cartridge dispenser can pop the seals on the cartridge and cause leakage. Avoid rapid pumping of the cartridge handle and take care to slowly express the material out of the mix tip.

DO NOT REMOVE THE RED TAB ON THE DUALINE DISPENSER. Removal can create pressure on the Dualine cartridge, which can pop the seals and cause leakage.

1. Insert a Dualine cartridge into the dispenser and express a small amount to verify that both catalyst and base pastes are flowing. Place the mixing tip on

the cartridge and continue to express the material to ensure complete mixing. Discard the first 1/2" of material.

- Express DualLine directly into the denture or saddle areas. This process takes about 30 seconds for a full upper denture. Alternatively, express the material into a disposable plastic cup; then pour into the denture or load with a spatula. NOTE: The pouring or spatula techniques will increase the risk of trapping air bubbles in the reline.

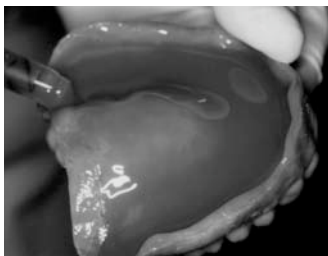- Wait about 2 - 3 minutes from the start of the dispensing before inserting the denture in the mouth. The time may vary slightly depending on ambient temperature. Wait 2 minutes initially; then adjust to your needs. Dry the tissue before insertion, if excessive saliva is present.
- Re-establish VDO, A-P and occlusal relationships while using a closed mouth impression technique and very light pressure. Excessive pressure will cause distortion in the tissue stop areas. Follow your normal border molding techniques. (The material has a flow similar to a light body

wash impression material.) Monitor the palate for excess material.

With removable partial dentures, insert and seat framework into the existing tooth rest areas to capture a nonfunctional impression. If a functional impression is desired, use TRIAD Regular Reline or TRIAD Hi-Flow Reline Material.

- Check the border areas for set before removing the denture at about 5 minutes from the start of dispensing DualLine. This will vary with ambient and intraoral temperatures.

Remove the relined denture carefully. Remember, the bond to the reline material is not complete until light cured in the TRIAD Unit.

Curing and Finishing

- Evaluate the reline and remove excess material from the denture. If voids or border shortages are apparent, hand-mix small amounts of DualLine for repairs or apply small amounts of TRIAD Regular Reline or TRIAD Hi-Flow Reline Sheet materials. Additional TRIAD Bonding Agent may be necessary at this time to facilitate bonding additional reline material.

If the impression is unacceptable, a lab knife may be used to remove the DualLine reline material since the bond is incomplete at this stage. The DualLine reline procedure can then be repeated.

- Carefully remove excess reline material with a scalpel, lab

knife or handpiece, as necessary, prior to light curing.

- Apply TRIAD Air Barrier Coating (ABC) to the DualLine surface and *cure in the TRIAD unit for 5 minutes*. Remove ABC completely with a soft brush or gauze sponge and warm soapy water. Evaluate fit.
- If a posterior seal is needed, slightly reduce the posterior palatal area by grinding. Apply and cure TRIAD Bonding Agent. Place TRIAD Regular Reline or TRIAD Hi-Flow Reline and shape to desired form. *Cure the posterior or seal in the TRIAD unit for 10 minutes*.

For a functional postdam, place in the mouth before curing and have the patient apply firm pressure for several seconds. Remove excess material and *cure in the TRIAD unit for 10 minutes*. Thin dentures, par-

tials with with metal frameworks and unknown base materials should be cured using several shorter cycles. Cool in water between cycles. Thermoplastic denture bases may be distorted in the TRIAD Unit.

- Finish, polish, and adjust occlusal relationships using conventional methods. Be sure to examine the interproximal areas for excess material. An explorer or scaler will remove excess easily, especially if a lubricant was used.

Laboratory Use

DualLine may be used as a laboratory reline material following standard clinical procedures to obtain an impression in the denture. Pour a cast and use the TRIAD Reline Fixture.

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